

Improve the quality of samples with a versatile medium

VIRCELL TRANSPORT MEDIUM



Transport Medium for the collection and preservation of **Viruses, Chlamydia and Mycoplasma**

MAIN BENEFITS

- Versatile medium, valid for different types of samples, microorganisms and diagnostic methods.
- Avoids the drying up of the sample.
- Improves microorganisms viability between collection and inoculation.
- Decreases bacterial overgrowth in contaminated samples.
- Successfully used for molecular based assays.
- Room temperature viability.
- Customized solution, available in different formats.
- Easy and safe to use, reduces the risk of contamination.



VIRCELL TRANSPORT MEDIUM


Improve the quality of samples
with a versatile medium

The success of virus diagnosis depends largely on the quality of the specimen and the conditions under which the specimen is transported and stored before being processed in the laboratory.

Vircell Transport Medium is intended for the collection and transport of clinical specimens containing viruses, chlamydia and mycoplasma from the collection site to the testing laboratory.

ASSAY PROCEDURE

Collect the sample with the sterile swab and put the swab into the tube containing the transport medium. Cut off the end of the swab stick standing out of the tube, taking care not to touch the rim of the tube. Close the tube tightly, and fill in the data of the patient, type of sample, and date and hour of collection. Send to the laboratory for analysis at room temperature. For prolonged transport, use a refrigerated delivery service at 2-8°C.

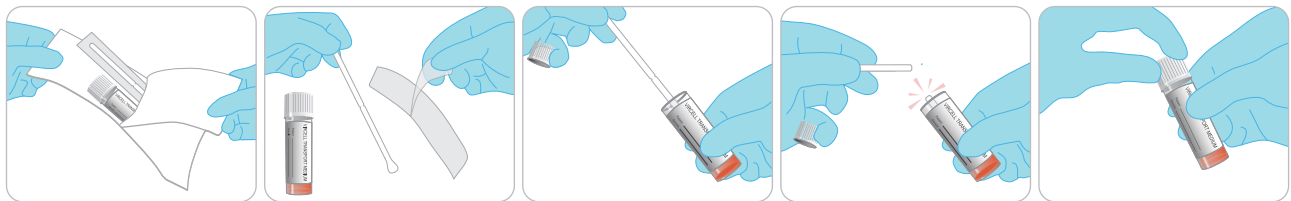

KIT CONTENTS

Vircell Transport Medium is provided in labeled screw-cap tubes - with different presentations - designed to preserve and transport the clinical sample. Some of our references incorporate 1 or 2 sterile swabs to ease specimen collection. Vircell medium consists of Hank's BSS containing HEPES, gelatine, BSA, sucrose and compatible antibiotics.

Cat. No.	Product description	Commercial swab	Pack size
TM011	TRANSPORT MEDIUM - 2 ml	1 sterile Rayon swab	50 units
TM003	TRANSPORT MEDIUM - 2 ml	1 sterile Dacron swab	50 units
TM005	TRANSPORT MEDIUM - 2 ml	No swab included	50 units
TM012	TRANSPORT MEDIUM - 2 ml	1 sterile nylon microfiber swab	50 units
TM004	TRANSPORT MEDIUM - 100 ml	No swab included	1 unit

PME010-05/20

vircell 
MICROBIOLOGISTS

Distributed by Abacus dx

1800 ABACUS (AUS) 0800 222 170 (NZ) | info@abacusdx.com | www.abacusdx.com

abacus dx